

## Micro cantilever

### Product name

OMCL-AC55TN-B3

Micro Cantilever	
<b>OmegaLevel</b> OMCL-AC55TN-B3	
Lot No.	
Nominal value	Inspection result
Resonant frequency 1.6 (MHz)	-
Spring constant 85 (N/m)	(Calculated value)
Inspector	
<b>OLYMPUS</b> Made in Japan	

Olympus Mega Hz Silicon cantilever with a sharpened tetrahedral tip


### OMCL - AC 55 T N - B3

OMCL : Olympus Micro Cantilever  
 AC : Main application is AC mode measurement  
 55 : Lever length of 55  $\mu\text{m}$   
 T : Sharpened Tetrahedral tip  
 N : No reflex metal coating  
 B3 : 18 chips / unit  
 Rectangular cross section chip

### Chip

The chip has a pentagonal cantilever on beveled chip shoulders side.

#### Dimensions

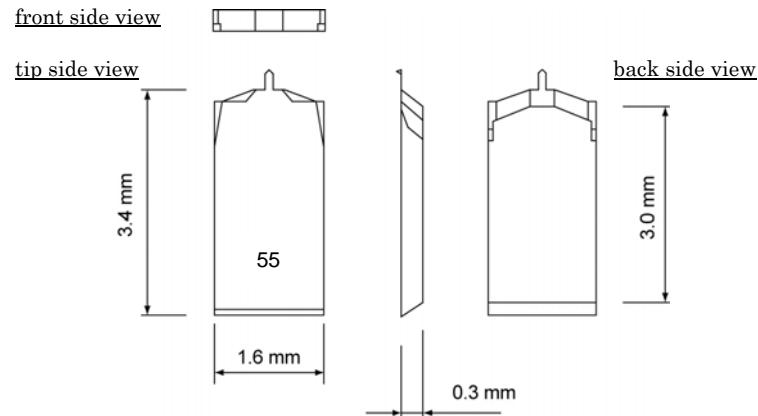


### Material

Tip & Lever	Silicon (n-type, 0.01 – 0.02 ohm.cm)
Metal coating (tip side)	Non
Metal coating (reflex side)	Non
Chip	Silicon (n-type, 0.01 – 0.02 ohm.cm)

### Probe

The probe is a sharpened tetrahedral. The probe is fabricated on the exact end of each cantilever.

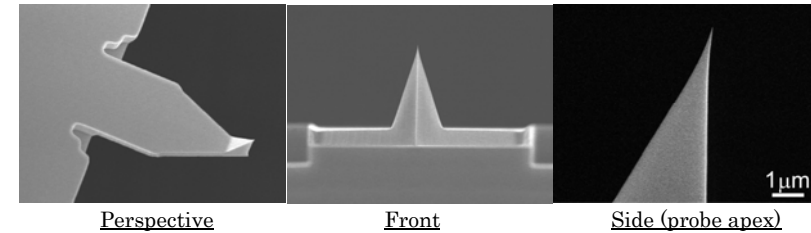


#### Dimensions

	Nominal value	Typical range
Probe length ( $\mu\text{m}$ )	12	8 - 16
Tip radius (nm)	7	4 - 10
Probe tip half angle (deg.)	(axis) less than 17.5 (side) less than 17.5	
Probe side tip angle (deg.)	(front) 0, (back) 35 (side) 18, 18	

### Cantilever

#### Dimensions

Cantilever length L ( $\mu\text{m}$ )	55 ( $\pm 10$ )
Cantilever width W ( $\mu\text{m}$ )	31 ( $\pm 1$ )
Cantilever thickness t ( $\mu\text{m}$ )	2.35 ( $\pm 0.3$ )

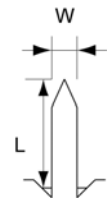

#### Calculated mechanical properties

	Nominal value	Typical range
Resonant frequency (MHz)	1.6	0.85 - 2.5
Spring constant (N/m)	85	38 - 184

**OLYMPUS**

OMCL-AC55TN-B3

Ver. 1.0 Jul 1, 2011